



96-Well CacoReady

User's Manual

LIMITED SINGLE-USE LICENSE

This is a binding legal agreement (this "License") between you (the "Customer" or "You"), and MedTech Barcelona regarding the enclosed "CacoReady™" System (this "Product"). Carefully read this License before opening the sealed package. BY OPENING THE SEALED PACKAGE, YOU AGREE TO BE BOUND BY THE TERMS OF THIS LICENSE. If You do not agree to the terms of this License, and You have not opened the sealed package, promptly notify MedTech Barcelona in writing and return the unopened Product and all accompanying items (including all written material and containers) to MedTech Barcelona within five (5) business days.

This License sets forth the Customer's rights to use this Product. The meaning of "Product" hereunder shall include any parts thereof and any documentation made available for use with this Product. Upon receipt of the applicable fee, MedTech Barcelona grants to the Customer the personal, non-exclusive, and non-transferable right to use this Product worldwide. The Customer is solely responsible for the use of this Product by any person.

The Customer acknowledges and agrees that the manufacture, use, sale, or import of this Product may be subject to one or more issued or pending patent applications and the corresponding foreign equivalents owned or controlled by MedTech Barcelona. The Customer acknowledges and agrees that the documentation and components included in the Product are the proprietary and copyrighted material of MedTech Barcelona. The purchase and use of this Product shall be governed under MedTech Barcelona's standard Terms and Conditions, and any other written agreement between the parties, and such terms and agreements are incorporated herein by reference. Your purchase of this Product, subject to this License, conveys to the Customer the non-transferable right to use the purchased amount of the Product and components of the Product in this shipment in research conducted by the Customer; such right shall expire upon completion of the experiments stemming from this shipment.

The Customer shall not detach living cells from the Product and shall not amplify genetic material (DNA, RNA) for purposes other than quantitative or qualitative analysis of proteins or nucleotide sequences.

The Customer shall not sell or otherwise transfer (a) this Product, (b) any components of this Product, or (c) any materials made using this Product or its components. Also, the Customer shall not use this Product or its components for Commercial Purposes other than those defined herein, and especially will not be able to: (1) make use of this Product in manufacturing; (2) make use of this Product for therapeutic, diagnostic or prophylactic purposes; or (3) resale this Product, whether or not such Product is resold for research purposes. For clarification, the Customer may use and transfer information or materials made through the use of this Product to support the Customer's research in any manner the Customer deems necessary or appropriate, and the Customer may provide services to third parties.

If you have any questions or need information on obtaining a license to use this Product or any of its components for purposes beyond the scope of this License, please contact MedTech Barcelona at:

reagents@medtechbcn.com

Table of contents

<i>PRODUCT DESCRIPTION</i>	4
<i>INTENDED USE</i>	4
<i>PRINCIPLE</i>	4
<i>TIMELINE FOR DELIVERY AND EXPERIMENTAL PROCEDURES</i>	5
<i>EQUIPMENT (NOT INCLUDED)</i>	6
<i>CONSUMABLES</i>	6
<i>SOLUTIONS (MAY BE INCLUDED)</i>	6
<i>HANDLING</i>	7
Replacement of Shipping Medium	7
<i>ASSAY PERFORMANCE</i>	7
<i>QUALITY CONTROLS</i>	8

Product Description

CacoReady is an *in vitro* cell-based model built on 21-day seeded human colorectal adenocarcinoma-derived cells (Caco-2), which resemble the intestinal epithelial barrier.

Caco-2 cells are considered a reference cellular model for assessing compound intestinal toxicity and their mechanisms of action. Furthermore, CacoReady can easily be adapted for high-throughput screening in drug discovery.

CacoReady is delivered in a 96-well/plate format with a unique Shipping Medium (a gel-like cell culture medium) established by MedTech Barcelona which enables the transport of cells at room temperature and in a ready-to-use format.

Intended Use

This product is mainly indicated for *in vitro* evaluation of:

- intestinal toxicity (e.g., cytotoxicity, genotoxicity, apoptosis, oxidative stress)
- mechanisms of action
- inflammatory and anti-inflammatory activities
- others

NOTE: This cell-based model is intended for scientific research purposes only. Not for human or veterinary use.

Principle

Caco-2 cells are an immortalized cell line derived from a human colorectal adenocarcinoma. These cells spontaneously differentiate into a monolayer that exhibits specific properties of absorptive intestinal epithelial cells, such as a well-defined brush border and intercellular junctions.

This cell line seeded on 96-well plastic plates is widely used to screen for intestinal adverse drug effects. Although limitations such as their tumor origin and the absence of the intestinal physiological context must be considered when extrapolating the data obtained *in vitro*, this cell line has the advantage of being easy to handle, reproducible, and adaptable to automatic high-throughput screening, which makes them relevant as an excellent toxicokinetic tool for predicting intestinal adverse drug effects and mechanism-based toxicity.

In the experimental setup, these cells are seeded for 21 days on 96-well plates (FALCON Cat# 353072) to form a tight cell monolayer.

Handling and experimental procedures are provided below. The manual has been written for users with experience in cell culturing and pharmacological drug discovery *in vitro* testing experiments. For more detailed advice and training opportunities, please contact us at reagents@medtechbcn.com.

Timeline for Delivery and Experimental Procedures

- Day 1: Start of Production (Seeding of cells)
- Days 14–15: Package Dispatch (depending on destination)
- Days 16–17: Package Delivery
- Day 18: Replacement of Shipping Medium (liquefaction)
- Day 21–25: Medium Replacement/Assay Performance (see Table 1)

Packages are dispatched on Mondays/Tuesdays and delivered within 24–48 h to EU countries, 48–72 h to USA, and 48–96 h to Asian countries. For other locations and customized schedules, please contact us at reagents@medtechbcn.com.

The recommended timing overview for permeability assays is Day 21 (Monday) (see Figure 1 for details).

CacoReady	Monday	Tuesday	Wednesday	Thursday	Friday
Week 0	12:00 p.m. (CET) last ordering day				
Week 1	Pre-Production		Start of Production Day 1		
Week 2	Day 7				
Week 3	Shipment Day 14	Reception of Plates			Liquefaction Day 18
Week 4	Perform Assay				Day 25
	Day 21				

Figure 1. Timeline of manufacturing and operation for CacoReady in 96-well format.

In case you choose to conduct the experiments later, the medium replacement must be carried out as follows:

DAY OF EXPERIMENT	MEDIUM REPLACEMENT
Tuesday (Day 22)	Monday (Day 21)
Wednesday (Day 23)	Monday (Day 21)
Thursday (Day 24)	Monday (Day 21) Wednesday (Day 23)
Friday (25)	Monday (Day 21) Wednesday (Day 23)

Table 1. Recommended day for medium replacement.

NOTE: These steps enable the planning of the assay according to the user's convenience.

Equipment (not included)

- Cell culture laminar flow hood
- CO2 incubator
- Water bath
- Multichannel pipettes
- Automatic multichannel micropipette (recommended)
- Aspiration system
- 96-well format vacuum manifold (Drummond Cat# 3-000-093 recommended)

Consumables

- Sterile culture medium containers (i.e., Costar 50 ml, Cat# 4870) (*not provided*)
- 15 and 50 mL conical tubes and 1.5 mL Eppendorf tubes (*not provided*)
- Pipette tips (*not provided*)

Solutions (may be included)

NOTE: MedTech Barcelona can supply Medium if required.

- Caco-2 Cell Culture Medium: Dulbecco's Modified Eagle's Medium – low glucose (1 g/L) (SIGMA Cat# D5546) supplemented with (final concentrations):
 - 10 % V/V Fetal Bovine Serum (BIOWEST Cat# DE14-801F)
 - 2 mM L-glutamine (LONZA Cat# BE17-605F)
 - 100 U/mL; 0.1 mg/mL Penicillin-Streptomycin (LONZA Cat# DE17-602F)

NOTE: If the specified reagents are not available, other reagents with similar features and specifications can be used.

Handling

Upon reception, retrieve the zipped bags containing the plates. Open the zip and leave the bag in a dark location at room temperature until Day 18 (refer to Timeline; Figure 1).

Replacement of Shipping Medium

CAUTION: *Never handle more than one plate at a time while changing the shipping medium. Re-solidification of the shipping medium may damage the cell monolayer.*

These steps will be carried out on Day 18 (refer to Timeline; Figure 1). Perform all manipulation under sterile conditions.

1. Retrieve the plates from the bags and remove the parafilm wrap.
2. Incubate the plates in a 5 % CO₂ humidified atmosphere at 37 °C for 4 hours, until the shipping medium reaches liquefaction.
3. Remove one CacoReady plate from the incubator and place it inside the laminar flow hood.
4. Using sterile procedures (inside the laminar flow hood), fill a sterile reagent reservoir with 20 mL of pre-warmed (37 °C) Caco-2 cell culture medium.
5. Open the CacoReady plate and leave the lid upwards, next to the plate.
6. Using the aspiration manifold connected to a vacuum pump (adjust aspiration flux to medium-low), aspirate the liquefied shipping medium, taking care not to disrupt the monolayer. Make sure the shipping medium has been removed from all wells. A remanent medium will be left in each well.
7. Using a multichannel pipette, dispense 150 µL of Caco-2 cell culture medium from the sterile reservoir, and fill the 96-wells of the CacoReady plate, column by column. Always add the medium against the wall of the well and not directly onto the cell monolayers.
8. Replace the lid and place the plate inside the cell culture incubator, set at 37 °C and 5 % CO₂.
9. Once the shipping medium has been substituted by fresh Caco-2 cell culture medium, the plates should be placed inside the incubator until next Monday (Day 21). Replacement with a new fresh medium will be carried out once (Day 21) or twice (Days 21 and 23), depending on the day of the assay (refer to Table 1 for details), following the procedure described above.

Assay Performance

Perform the assay as needed. Please, feel free to contact us with any queries or doubts.

Quality Controls

This section provides the morphology of 21-day differentiated Caco-2 cells seeded on plastic dishes and the pre-specified criteria of monolayer confluency for batch release and after recovery from the shipping medium (SM).

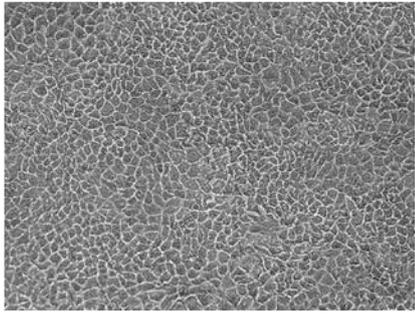
Cell morphology	Monolayer confluency	
	Batch release	Recovery from SM
	$\geq 95 \%$	$\geq 90 \%$

Table 2. Morphology of 21-day differentiated Caco-2 cells.